

# JUNIOR

Biologics Formulation



**UN**CHAINED  
LABS

## End old routines

Junior automates all your hands-on, back and forth and hard to do biologics formulation tasks. It ties them together and makes them one totally routine, high-throughput process. Hand off the grunt work to Junior, and hit the next big thing on your plate instead.

**Characterize formulations**

**Developability assessments**

**Evaluate formulation robustness**

**Preformulation and formulation development**



## Hand off your process

Walk up, set up your run and walk away. Junior automates the manual measurements you do one by one at the bench. Choose a combo of high-throughput pH with visual inspection, visual inspection with viscosity or viscosity with pH. Crank through more samples in a single day, get them done the same way every time and check out a broader developability and formulation space.



- |   |                     |   |                           |
|---|---------------------|---|---------------------------|
| 1 | Vial/plate gripper  | 5 | 3-Position plate rack     |
| 2 | Waste bin           | 6 | Visual inspection station |
| 3 | Viscosity station   | 7 | Vial/plate hotel          |
| 4 | Tool and plate rack |   |                           |

## Bench the pH meter


One at a time pH measurements just don't cut it. A high-resolution, 4-probe pH arm automates the whole process. It blasts through a 96-well plate in less than 30 minutes – that includes calibrations, measurements, recording and washing. So no more prep, no more mess, and no more standing at the bench.



Multi-channel pH probe

## Knock viscosity out in 6


Get viscosity measurements for Newtonian and non-Newtonian fluids in under 6 minutes with as little as 100  $\mu\text{L}$  of sample. With its spot-on temp control, the Junior viscosity station gives you the consistency and accuracy you need. Automated washing kills any carryover between samples. Crank through your viscosity measurements 5X faster, and cut your hands-on time down by 40X.

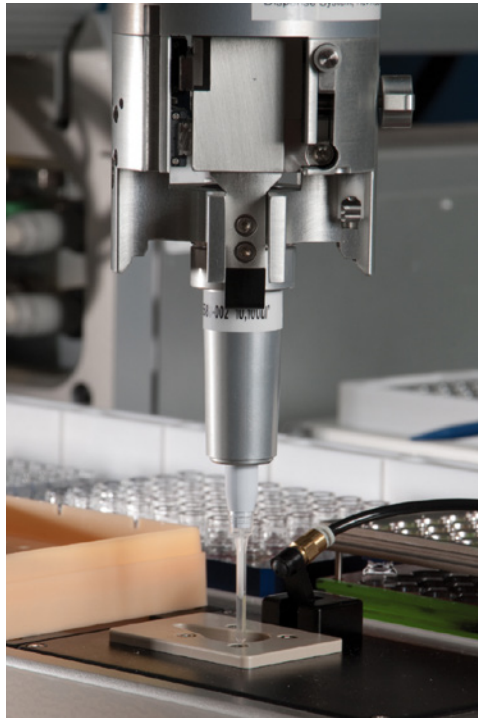

Viscosity station

## Dead-on visual inspection


Junior does a 3-in-1 visual inspection using a fast, non-destructive assay. Grab color, turbidity and count visible particles all in one shot. It also takes a snapshot of every sample and matches those up against standards. That means you get a more accurate analysis without all the subjectivity. Images are permanently archived so you can go back and check them out anytime.

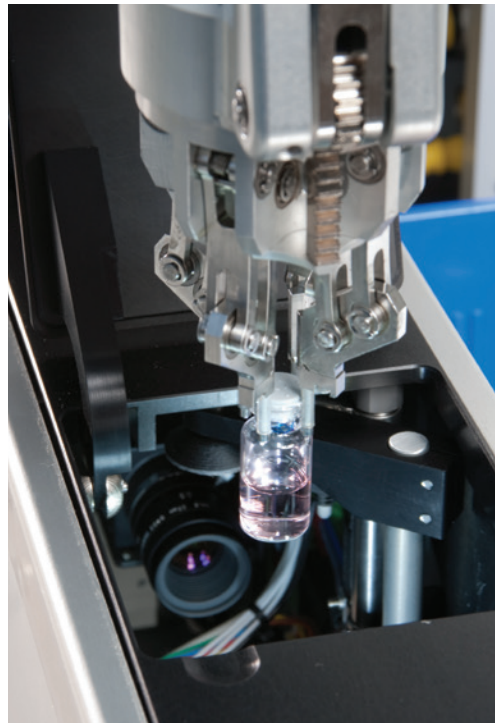

Visual inspection station

## Break through bottlenecks

LEA doesn't move the bottleneck, it totally unclogs it. You get intuitive experimental design that makes sense on the front end and full integration with analytical tools on the back end. LEA also links your conditions, steps and analytical data together so your final report doesn't just have numbers, it's got all the information you need to make real decisions.


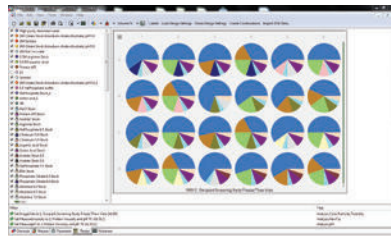


### Library Studio

Design complex, high-throughput experiments in an array-based format

### PolyView

Review and report all info from experimental design, execution and analytics

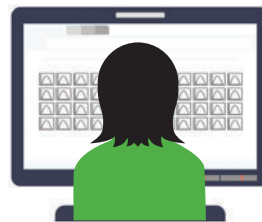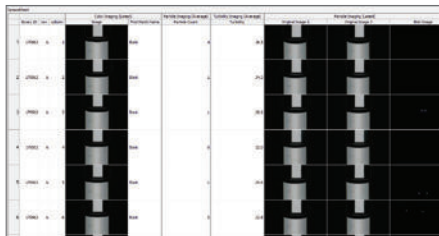

### Automation Studio

Execute designed experiments and integrated analytics



**Unchained Labs**

6870 Koll Center Parkway

Pleasanton, CA 94566

Phone: 1.925.587.9800

Toll-free: 1.800.815.6384

Email: [info@unchainedlabs.com](mailto:info@unchainedlabs.com)

© 2017–2018 Unchained Labs. All rights reserved. Junior is a trademark and Unchained Labs is a registered trademark of Unchained Labs. All other brands or product names mentioned are trademarks owned by their respective organizations.